




Bürkle® Discharge hopper single, V2A

Product Images

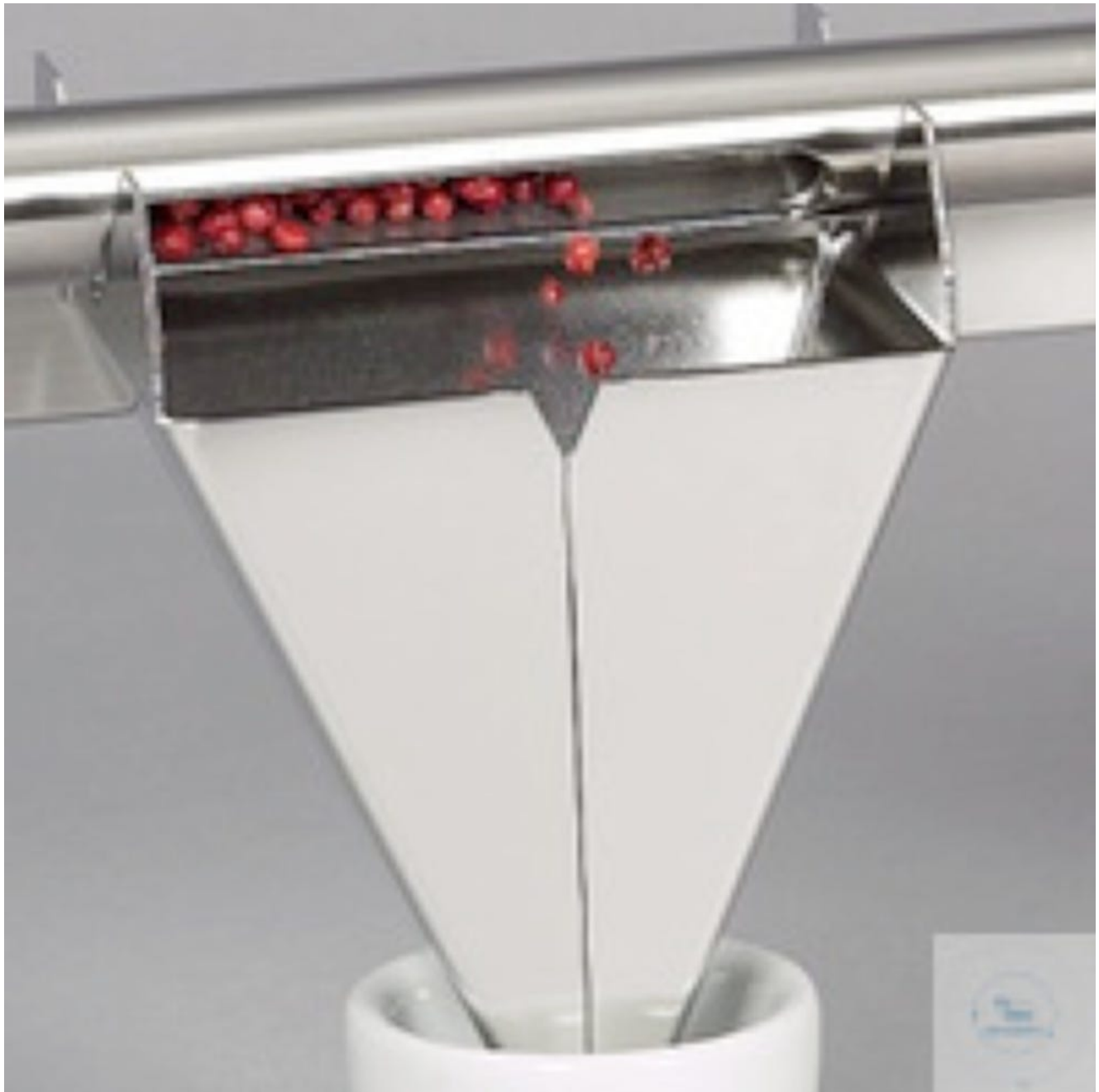

Description

The separate emptying of the multi-collector chambers can be done by means of Quali-Paper or closure sleeves. It is even more convenient and also much faster with the emptying rack. The samples in the collection chambers can be emptied directly into the sample vessels. Simply place the multicollector on the emptying rack, place the sample containers under the funnels, open the chambers - and you're done. This saves valuable time, nothing contaminates the samples, they are not mixed and nothing is lost. The individual funnels can be variably adjusted to the different center distances of the collection chambers. Thus, the emptying frame fits the multi-collectors with 550 mm and 850 mm length. If several emptying racks are bolted together, multicollectors of any length can be emptied. Emptying device, complete with three hoppers Stainless steel V2A (1.4301) 100 cm long

Additional Information

No.	KL-0715
Manufacturer (Brand)	Bürkle
VGKL number	100450033

