




neolab Migge GmbH
Rischerstr. 7-9
69123 Heidelberg
Deutschland
+49 (0)6221 /
8442-44
<https://www.neolab.de>
e

Umsatzsteuer-
Identifikationsnummer
:
DE 143 450 657

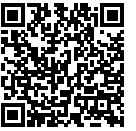

SERVA BlueClear SB Stripping buffer 250 ml

€47.00
plus VAT &
Shipping

Product Images



Description

Ready-to-use buffer for efficient stripping of high-affinity antibodies from Western blot membranes after chemiluminescence detection.

Easy to use - incubating the membrane in the stripping buffer for 30-60 min at room temperature is sufficient.

Hard to remove antibodies are efficiently stripped by incubation in heated buffer.

After washing in PBST or TBST, the membrane can be blocked again and incubated with antibodies.

- Fast, simple and highly efficient
- Gentle
- Without β -mercaptoethanol or DTT
- Suitable for nitrocellulose and PVDF membranes

pH: 2.6-2.8

Additional Information

No.	SV-0096
Manufacturer (Brand)	SERVA
EAN	4058072372330
Transport temperature	Room temperature

