




Bandelin K 6, INSERT BASKET

Product Images

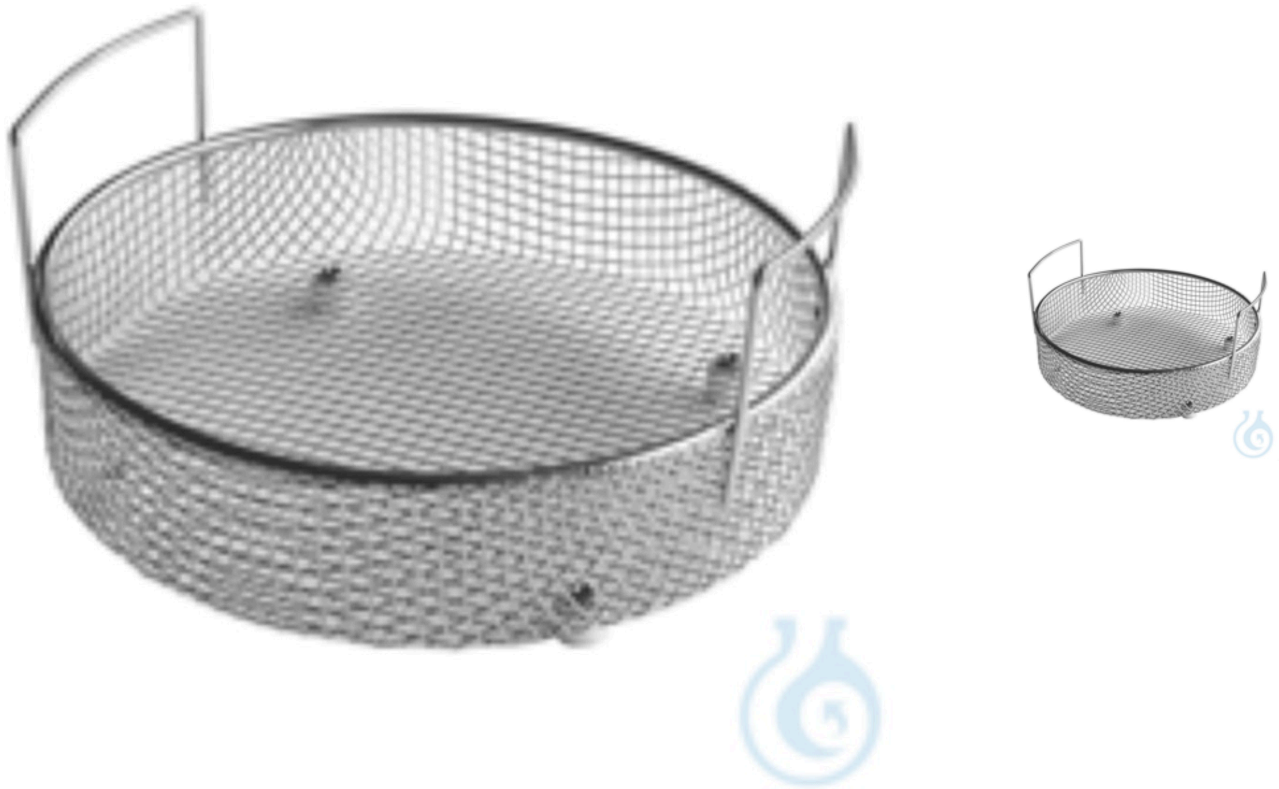

Description

For holding sonication material; made of stainless steel During ultrasonic operation, nothing may be placed on the bottom of the tank, as abrasion may occur between the sonication material and the ultrasonic oscillating tank. It is therefore recommended to use a suitable basket to hold the sonication material. The BANDELIN baskets are designed in such a way that the sonication material is always in the optimum range of the cavitation field (1-2 cm from the floor) in the oscillating tank. With the aid of a basket, the sonication material can also be conveniently placed in the unit and removed after operation, e.g. for rinsing, without the user coming into contact with the sonication liquid in the oscillating tank.

Additional Information

No.	BL-0220
Dual use indicator	false
Manufacturer (Brand)	Bandelin
VGKL number	960731103
EAN	4058072505806
Country of origin	Deutschland
Region of origin	Berlin
Packing width	0.22000 m
Packing height	0.07000 m
Packing depth	0.22000 m
Packaging volume	0.00332700 ccm
Customs tariff number	84799070

