




Bandelin SD 06, GLASBECHER KPL.

Product Images

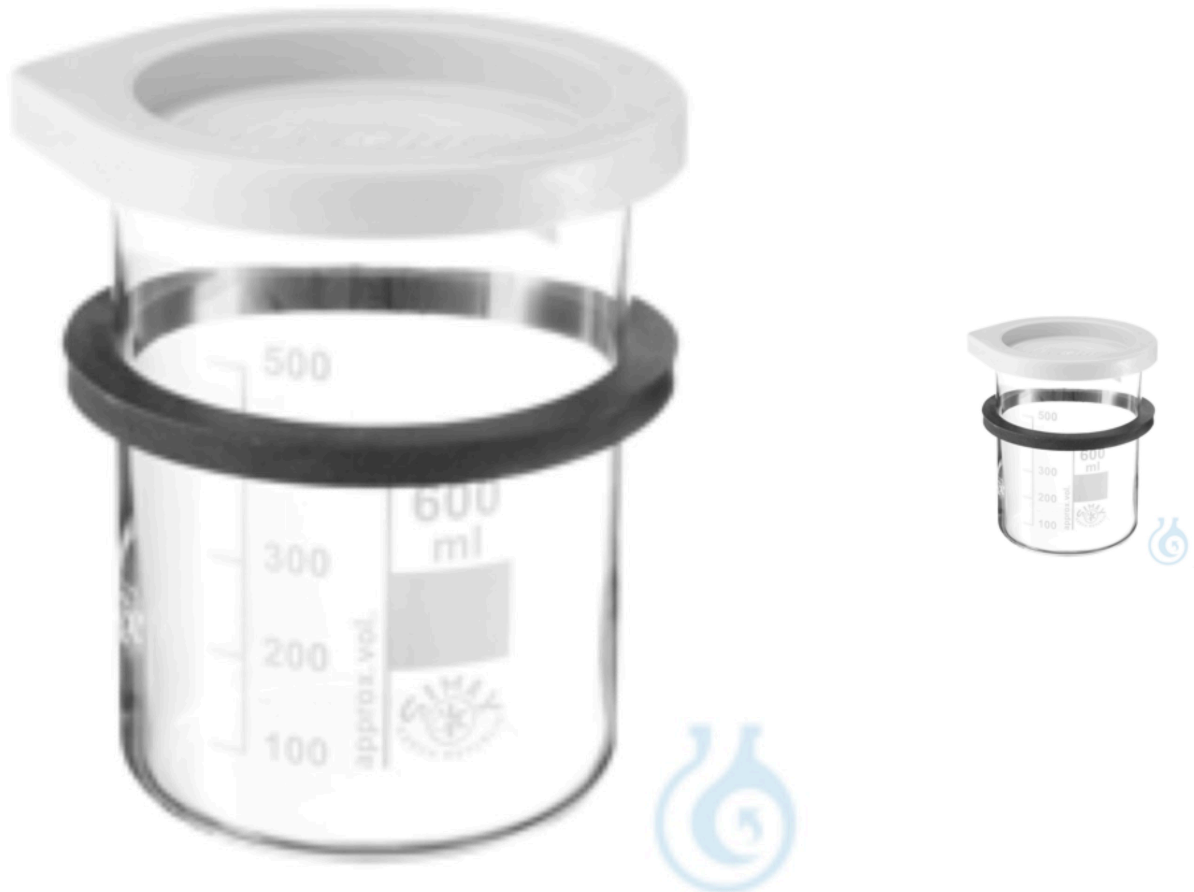

Description

Glass, 600 ml, for insertion into perforated lid/insertion strip with ring GR 06 and lid DD 06 Insertion beakers are used when applying chemically aggressive liquids (e.g. TICKOPUR J 80 U), removing chemically aggressive or abrasive contaminants or sonicating sample liquids (indirect sonication). Because nothing may be placed on the bottom of the ultrasonic bath, BANDELIN offers suitable perforated lids/insert strips to hold the insert beakers. Further advantage: Several insert cups can be sonicated with different liquids at the same time (e.g. removal of dental cements with STAMMOPUR Z in one and cleaning with simultaneous chemical disinfection of dental instruments with STAMMOPUR DR 8 in a second insert cup) Please note that for the transmission of the ultrasound to the insert cup the surface tension of the water must be reduced and therefore a cleaning solution (e.g. TICKOPUR R 33 - 3%) must be used as coupling liquid / contact liquid.

Additional Information

No.	BL-0193
Dual use indicator	false
Manufacturer (Brand)	Bandelin
VGKL number	960736101
EAN	4058072505387
Country of origin	Tschechien
Volume	600 ml
Packing width	0.11000 m
Packing height	0.18000 m
Packing depth	0.09000 m
Packaging volume	0.00119700 ccm
Wide	90 mm
Customs tariff number	70200080
Height	125 mm
Length	90 mm

